**EN Mediblink Blood Pressure Monitor Experience M550 AFib**

**INSTRUCTIONS FOR USE**

Please read the instructions for use carefully before using the product.

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# TABLE OF CONTENTS

|                                                                         |           |
|-------------------------------------------------------------------------|-----------|
| <b>Safety Notice</b>                                                    | <b>06</b> |
| Important Instructions Before Use                                       | 09        |
| <b>Device Illustration</b>                                              | <b>11</b> |
| Monitor Device                                                          | 11        |
| Display                                                                 | 12        |
| Content                                                                 | 13        |
| <b>Important Testing Guidelines</b>                                     | <b>14</b> |
| <b>Quick Start</b>                                                      | <b>15</b> |
| Lithium Battery Charging                                                | 16        |
| System Settings                                                         | 17        |
| Applying the Rigid Arm Cuff                                             | 19        |
| <b>Device Operation</b>                                                 | <b>20</b> |
| Testing                                                                 | 20        |
| Atrial fibrillation (AFib)                                              | 22        |
| Information for the doctor on frequent appearance of the AFib indicator | 22        |
| Power Off                                                               | 22        |
| Last 3 Tests Average                                                    | 23        |
| Memory Check                                                            | 23        |
| Memory Deletion                                                         | 24        |
| Low Battery Indicator                                                   | 24        |
| For Service Personnel – Static Pressure Measurement                     | 25        |
| Arm Movement Indicator                                                  | 25        |
| Loose Cuff Indicator                                                    | 25        |
| Troubleshooting                                                         | 26        |
| <b>Blood Pressure Information</b>                                       | <b>28</b> |
| Blood Pressure                                                          | 28        |
| Operating principles                                                    | 28        |
| WHO Blood Pressure Classification Indicator                             | 29        |
| Health Reminder                                                         | 30        |
| <b>Blood Pressure Q&amp;A</b>                                           | <b>31</b> |
| <b>Maintenance</b>                                                      | <b>32</b> |
| <b>Specifications</b>                                                   | <b>34</b> |
| <b>Electromagnetic Compatibility Information</b>                        | <b>36</b> |
| <b>Additional Notes</b>                                                 | <b>39</b> |
| <b>Warranty</b>                                                         | <b>41</b> |

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# Safety Notice

---

Thank you for purchasing the Mediblink Blood Pressure Monitor Experience M550 AFib. The device has been constructed using reliable circuitry and durable materials. Used properly, this device will provide years of satisfactory use.

## Target group

It is designed for self-measurement by adults in the home environment and is suitable for users whose upper arm circumference is within the range printed on the cuff. The blood pressure part is for adults and adolescents over 12 years of age, and the AFib part is for adults aged 18 and above. This was successfully tested as part of a clinical study.

## Indications for use

The device is designed for home or clinical use. All functions can be used safely and values can be read out in one LCD display. Measurement position is on adult upper arm only. Please read this manual thoroughly before using the device. Please retain this manual for future reference. For specific information about your blood pressure, please consult your doctor.

To avoid risk and damage, follow all warning precautions. Operate device only as intended. Read all instructions prior to use.

## Warning signs and symbols used

|                                                                                       |                                                                                           |                                                                                   |                                 |
|---------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|---------------------------------|
|      | Caution                                                                                   |  | Serial number                   |
|      | Mandatory                                                                                 |  | Medical device                  |
|      | Prohibited                                                                                |  | Temperature limitation          |
|      | Type BF equipment                                                                         |  | Humidity limitation             |
|      | Refer to instruction manual                                                               |  | Atmospheric pressure limitation |
| <a href="http://www.mediblink.com/t/m550.pdf">http://www.mediblink.com/t/m550.pdf</a> |                                                                                           |  | Reference number                |
|       | Authorized Representative in European Union                                               |  | Model number                    |
|      | Number of products included in packaging                                                  |  | Batch number*                   |
|      | General warning                                                                           |  | Unique Device Identifier        |
|       | CE Marking of Conformity, Notified body ID number                                         |  | Recyclable                      |
|      | Discard the used product to the recycling collection point according to local regulations |  | Paper                           |
|      | Manufacturer                                                                              |  | Keep dry                        |
|      | Date of manufacture                                                                       |  | Keep away from sunlight         |

\*The production date can be read from LOT number [YYYYMM]; first four digits represent the Year and last two digits Month of production. Example: LOT 202503 = March 2025

### Caution

Individuals with serious circulation problems may experience discomfort. Consult your physician prior to use.

Contact your physician if test results regularly indicate abnormal readings. Do not attempt to self-treat these symptoms without consulting your physician first.

The product is designed for its intended use only. Do not misuse in any way.

The product is not intended for infants or individuals who cannot express their intentions.

Do not disassemble or attempt to repair.

Do not use cell phones and other devices, which generate strong electrical or electromagnetic fields, near the device, as they may cause incorrect readings and interference or become a source of interference to the device.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the Competent Authority of the Member State in which the user and / or patient is established.

### Battery Precautions

This product contains lithium battery, do not discard as normal waste.

Recharge batteries when Low Battery Indicator “” appears on screen.

Recharge the lithium battery at regular intervals to ensure proper device functionality.

Dispose batteries properly; observe local laws and regulations.

Users are not allowed to replace the lithium battery by themselves. If it is confirmed that the lithium battery is faulty, please do not repair it by yourself. Please contact the Manufacturer or retailer.

## Important Instructions Before Use

1. Do not confuse self-monitoring with self-diagnosis. Blood pressure measurements should only be interpreted by a health professional who is familiar with your medical history.
2. Contact your physician if test results regularly indicate abnormal readings.
3. If you are taking medication, consult with your physician to determine the most appropriate time to measure your blood pressure. NEVER change a prescribed medication without first consulting with your physician.
4. Individuals with serious circulation problems may experience discomfort. Consult your physician prior to use.
5. For persons with irregular or unstable circulation resulting from diabetes, liver disease, arteriosclerosis or other medical conditions, there may be variations in blood pressure values measured at the wrist versus at the upper arm. Monitoring the trends in your blood pressure taken at either the arm or the wrist is nevertheless useful and important.
6. People suffering from vascular constriction, liver disorders or diabetes, people with cardiac pacemakers or a weak pulse, and women who are pregnant should consult their physician before measuring their blood pressure themselves. Different values may be obtained due to their condition.
7. People suffering from arrhythmias such as atrial or ventricular premature beats or atrial fibrillation should only use this blood pressure monitor in consultation with your doctor. In certain cases, oscillometric measurement method can produce incorrect readings.
8. Too frequent measurements can cause injury to the patient due to blood flow interference.
9. The cuff should not be applied over a wound as this can cause further injury.
10. **DO NOT** attach the cuff to a limb being used for IV infusions or any other intravascular access, therapy or an arterio-venous (A-V) shunt. The cuff inflation can temporarily block blood flow, potentially causing harm to the patient.
11. The cuff should not be placed on the arm on the side of a mastectomy. In the case of a double mastectomy use the side of the least dominant arm.
12. Pressurization of the cuff can temporarily cause loss of function of simultaneously used monitoring equipment on the same limb.
13. A compressed or kinked connection hose may cause continuous cuff pressure resulting in blood flow interference and potentially harmful injury to the patient.
14. Check that operation of the device does not result in prolonged impairment of the circulation of the patient.


15. The product is designed for its intended use only. Do not misuse in any way.
16. The product is not intended for infants or individuals who cannot express their intentions.
17. Prolonged over-inflation of the bladder may cause ecchymoma of your arm.
18. Do not disassemble the device or arm cuff. Do not attempt to repair.
19. Use only the approved arm cuff for this device. Use of other arm cuffs may result in incorrect measurement results.
20. The system might produce incorrect readings if stored or used outside the manufacturer's specified temperature and humidity ranges.
21. Do not use the device near strong electrical or electromagnetic fields generated by cell phones or other devices, they may cause incorrect readings and interference or become interference source to the device. Do not use the device during patient transport outside healthcare facility for interference source existing as well.
22. This product contains lithium battery, do not discard.
23. Recharge the lithium batteries when Low Battery Indicator "☐" appears on screen.
24. Recharge the lithium battery at regular intervals to ensure proper device functionality.
25. Dispose batteries properly; observe local laws and regulations.
26. Users are not allowed to replace the lithium battery by themselves. If it is confirmed that the lithium battery is faulty, please do not repair it by yourself. Please contact the Manufacturer or dealer.
27. Only use a recommended AC adaptor double-insulated complying with EN 60601-1 and EN 60601-1-2 (not included). An unauthorized adapter may cause fire and electric shock.
28. Advising operator that Instruction manual/Booklet must be consulted.
29. Essential performance:

|                                                            |                    |                |
|------------------------------------------------------------|--------------------|----------------|
| <b>Electrosurgery interference recovery</b>                | Refer 202.6.2.101  | IEC 80601-2-30 |
| <b>Limits of the error of the manometer</b>                | Refer 202.12.1.102 | IEC 80601-2-30 |
| <b>Reproducibility of the BLOOD PRESSURE DETERMINATION</b> | Refer 201.12.1.107 | IEC 80601-2-30 |


# Device Illustration

---

## Monitor Device

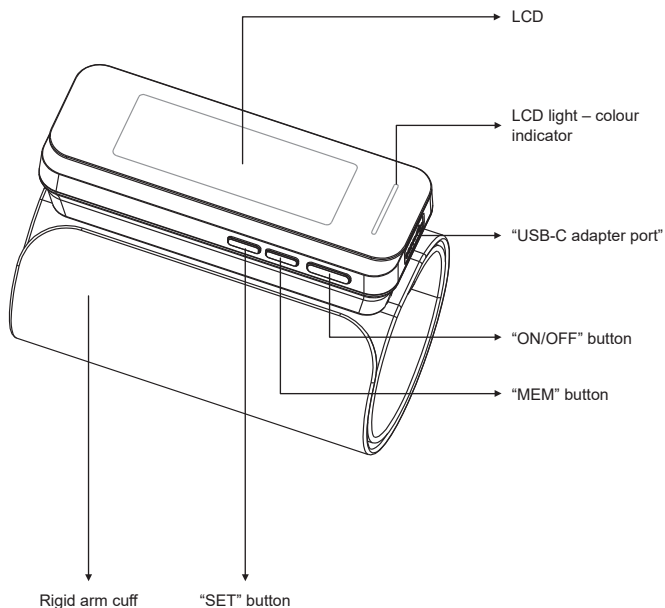
## Display
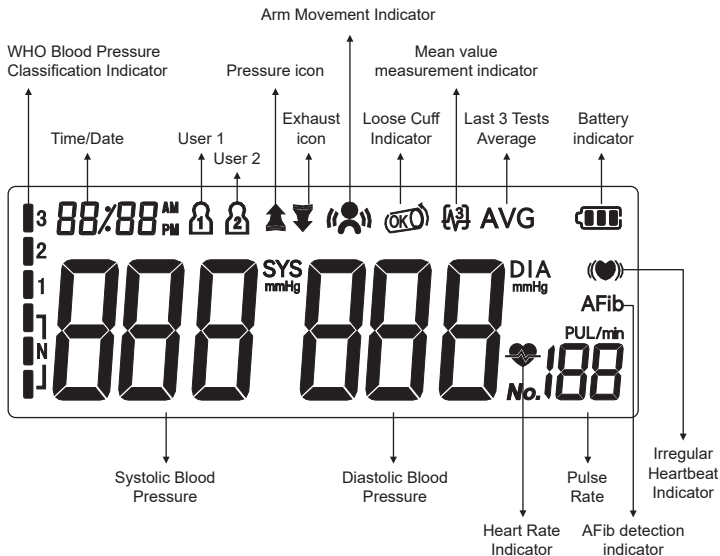
# Content


## 1. Blood Pressure Monitor Device

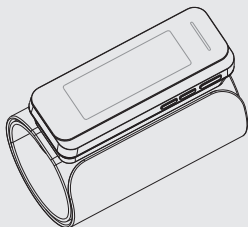
## 2. User's Manual

## 3. Storage Bag
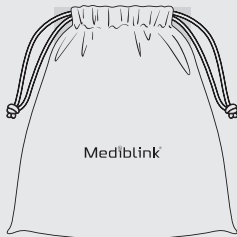
## 4. USB-C Cable
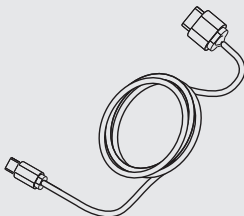
## Important Testing Guidelines

---

1. Avoid eating, exercising, and bathing for 30 minutes prior to testing.
2. Sit in a calm environment for at least 5 minutes prior to testing.
3. Do not stand while testing. Sit in a relaxed position while keeping your arm level with your heart.
4. Avoid speaking or moving body parts while testing.
5. While testing, avoid strong electromagnetic interference such as microwave ovens and cell phones.
6. Wait 3 minutes or longer before re-testing.
7. Try to measure your blood pressure at the same time each day for consistency.
8. Test comparisons should only be made when monitor is used on the same arm, in the same position, and at the same time of day.

- This blood pressure monitor is not recommended for people with severe arrhythmia.
- Do not use this blood pressure monitor if the device is damaged.

## Quick Start

---

1. Remove thick clothing from the arm area.
2. Rest for several minutes prior to testing. Sit down in a quiet place, preferably at a desk or table, with your arm resting on a firm surface and your feet flat on the floor. (See Figure A)

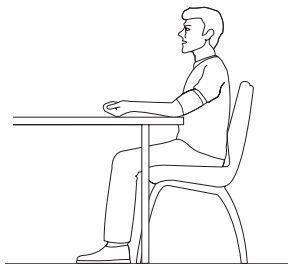
Figure A

3. Apply the cuff to your left arm and keep level with your heart. Bottom of cuff should be placed approximately 1–2cm (0.4–0.8") above elbow joint. (See Figures B and C)
4. Press "START/STOP" Button to start testing.

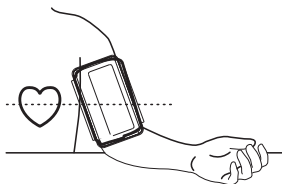
Figure B


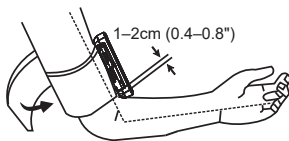
Figure C

## Lithium Battery Charging

The device has a built-in lithium battery that must be recharged using a specific power adapter. When low battery indicator icon “” appears on the screen, please connect the power adapter to charge. When charging, the screen will be displayed in a cycle from “” to “”. To charge the device, insert the adapter connector into the adapter interface (USB-A), then plug the adapter into the mains socket. Connect the charging cable to the USB-C adapter port on the device. Once connected to the power supply, the device will start recharging and the screen power display will change, indicating that you are using the adapter for charging or for power. The device can also be used to perform the measurement as usual by pressing the “ON/OFF” button while connected to power supply.

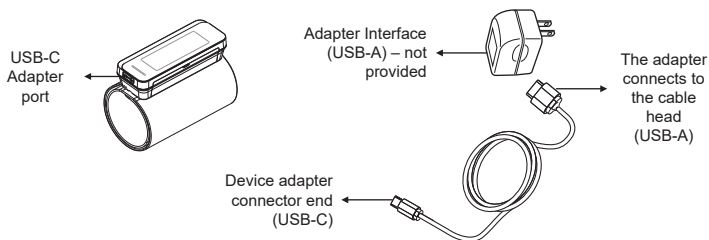The AC adapter port is located on the side of the monitor. A medical-grade 2 MOPP AC adapter (output: DC 5.0V, 1000mA) may be used with this device (recommended, not included). The DC plug must be a 2.1mm coaxial connector with center positive polarity (positive inside, negative outside). Do not use any other type of AC adapter, as it may damage the device. When connected to an external power adapter, the adapter is considered part of the medical electrical system.

**Note:** Power supply is specified as part of ME EQUIPMENT. Do not use other non-standard adapters. The above adapter parameters are as follows:

Input range: 100V–240V to 50/60Hz, 0.25A

Output range: 5.0V DC 1.0A, 5.0W


The interface at the adapter is USB-A port, and the interface at the device is USB-C port.

## System Settings

With power off, press “SET” button to activate System Settings. The Memory User icon flashes.

### 1. Select Memory User


While in the System Setting mode, you may accumulate test results into 2 different Users. This allows multiple users to save individual test results (up to 150 memories per User.) Press “MEM” button to choose a User setting. Test results will be automatically stored in each selected User.
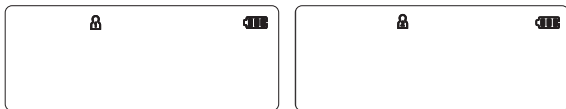
### 2. Time/Date setting


Press “SET” button again to set the Time/Date mode. Set the year first by adjusting the “MEM” button.

Press “SET” button again to confirm current month. Continue setting the date, hour and minute in the same way. Every time the “SET” button is pressed, it will lock in your selection and continue in succession (month, day, hour, minute, 12/24 hours).

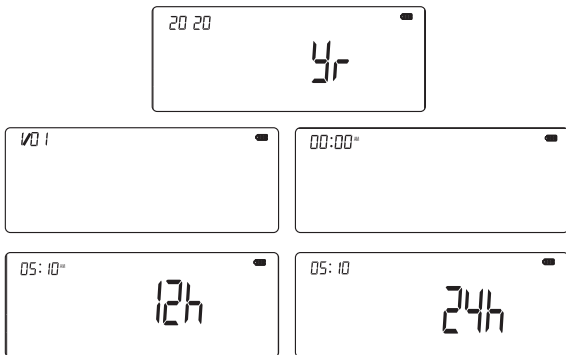
### 3. Time Format setting

Press “SET” button again to set the time format setting mode. Set the time format by adjusting the “MEM” button. EU means European Time, US means US Time.
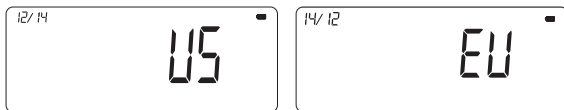
### 4. AFib function Setting

Press “SET” button to enter AFib function setting mode. Set AFib format ON or OFF by pressing the “MEM” button.
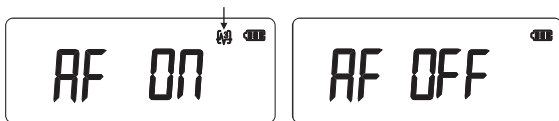
### 5. Mean value measurement mode time setting (only when AFib function is ON)

Press “SET” button to enter Mean value measurement mode time setting. Set Mean value measurement time format 15 seconds or 30 seconds or 60 seconds or 90 seconds by pressing the “MEM” button.


### 6. Saved Settings

While in any setting mode, press “ON/OFF” button to turn the device off. All information will be saved.



**Note:** If device is left on and not in use for 3 minutes, it will automatically save all information and shut off.

## Applying the Rigid Arm Cuff

1. Place the Blood Pressure monitor with open Rigid cuff on the upper arm and secure it firmly, ensuring proper fit and preventing air leakage.

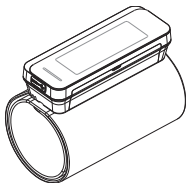2. With sticky nylon section facing outward, insert the end of the cuff underneath metal ring of cuff.
3. Fasten cuff about 1–2cm (0.4–0.8") above the elbow joint. For best results apply cuff to bare arm and keep level with heart while testing.

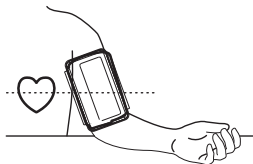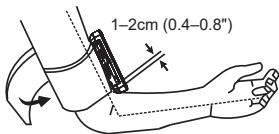
# Device Operation

## Testing

### 1. Power On


Press “START/STOP” button to turn the device on. The LCD screen will illuminate for one second as device performs a quick diagnosis. A long tone will indicate when device is ready for testing.

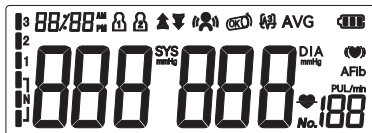**Note:** Device will not function if residual air from previous testing is present in cuff. The LCD will flash “↓” until pressure is stabilized.

### 2. Testing


After cuff inflation, air will slowly rise as indicated by the corresponding cuff pressure value. A flashing “♥” will appear simultaneously on the screen signaling heart beat detection. If the third consecutive measurement function is enabled (Mean value measurement mode), the first measurement starts and the measurement screen is displayed as shown in Figure b.

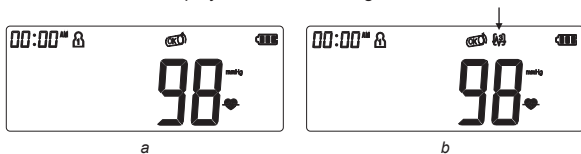**Note:** Keep relaxed during testing. Avoid speaking or moving body parts.


### 3. Result Display

- (1) Mean value measurement mode is **NOT** enabled.  
Four short beeps sound when testing is complete. The screen will display measurements for systolic and diastolic blood pressure. An indicator representing the current measurement will appear next to the corresponding WHO Classification.

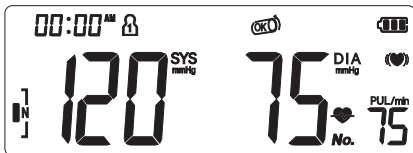(2) Mean value measurement mode is **ENABLED** (AFib function is ON).

When the Mean Value Measurement Mode is enabled, the device automatically performs three consecutive measurements. After each normal measurement, the result is displayed and a 15-second countdown appears on the screen before the next measurement begins. If a measurement error occurs, the corresponding error message is displayed and a 15-second countdown is performed before the next attempt. If all three measurements are unsuccessful, the sphygmomanometer will switch off automatically. If three valid measurements are completed, the device does not display the third individual result separately, but instead calculates and directly displays the average value of the three measurements. The indicator corresponding to the WHO blood pressure classification for the averaged result will appear on the display.

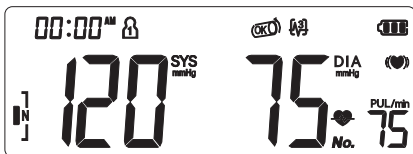**Note:** Refer to Page 29 for detailed WHO Blood Pressure Classification Information.

### Irregular Heartbeat Indicator

If the monitor detects an irregular heart rhythm two or more times during the measuring process, the Irregular Heartbeat Symbol “(♥)” appears on screen along with measurement results. Irregular heartbeat rhythm is defined as rhythm that is either 25% slower or faster than the average rhythm detected while measuring systolic blood pressure and diastolic blood pressure. Consult your physician if the Irregular Heartbeat Symbol “(♥)” frequently appears with your test results.

## Atrial fibrillation (AFib)

This device is able to detect atrial fibrillation (AFib). This symbol “**AFib**” indicates that atrial fibrillation was detected during the measurement. Please refer to the next paragraph for information regarding the consultation with your doctor.

### Information for the doctor on frequent appearance of the AFib indicator

This device is an oscillometric blood pressure monitor that also analyses pulse irregularity during measurement. The device is clinically tested. The AFib symbol is displayed after the measurement, if atrial fibrillation occurred during measuring. If the AFib symbol appears after having performed a full blood pressure measurement episode (triplicate measurements), the patient is advised to perform another measurement episode (triplicate measurements). If the AFib symbol appears again, we recommend the patient to seek medical advice.

If the AFib-symbol appears on the screen of the blood pressure monitor, it indicates the possible presence of atrial fibrillation. The atrial fibrillation diagnosis, however, must be made by a cardiologist based on ECG interpretation.


## Power Off

The “START/STOP” button can be pressed to turn off the device in any mode. The device will automatically switch off after approximately 3 minutes if no operation is performed in any mode.

**Safety Precaution:** If pressure in arm cuff becomes too extreme while testing, press the “START/STOP” button to turn power off. The cuff pressure will rapidly dissipate once the device is off.


## Last 3 Tests Average

With power off, press the “MEM” button to activate screen display. After the device performs a self-diagnosis, the screen will display the average test results from the last 3 readings of the last User used. The “AVG” symbol will appear along with the corresponding WHO Blood Pressure Indicator. The Memory Check mode can be accessed by pressing “MEM” button. To check the average results from other Users, select the desired User first prior to activating “SET” button in the off position (See “Select Memory User” on Page 17).

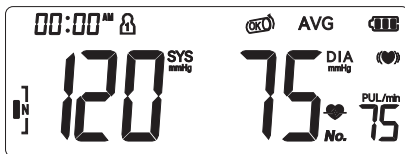
## Memory Check

You may check past test results by using the “MEM” button. The most recent test result and oldest test result in memory can be viewed by pressing and holding the “MEM” button. Upon activating test results, you can press the “MEM” button to scroll through all test results stored in memory.
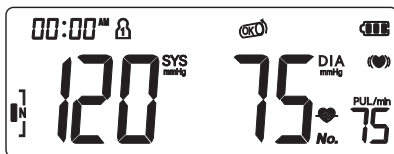
## Memory Deletion

Memory for a selected User may be deleted while in Memory Check mode. Press and hold the “SET” button for approximately 3 seconds to delete all memory records from the selected User. The monitor will beep indicating successful deletion and then transfer into testing mode. Press the “START/STOP” button to turn the device off.

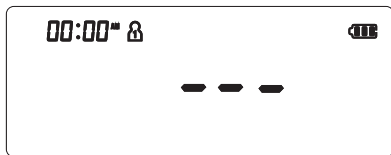**Note:** Memory cannot be recovered once it has been deleted.


## Low Battery Indicator

Four short warning beeps will sound when the battery power is low and insufficient to inflate the cuff for measurement. The “” symbol will appear on the display simultaneously for approximately 5 seconds before the device switches off automatically. Recharge the batteries at this time. No stored measurement data will be lost during this process.

## For Service Personnel – Static Pressure Measurement

In the power down state, press and hold the “START/STOP” button and until the LCD screen is fully displayed, then release the button. When the LCD screen displays the double zero, the blood pressure monitor is in static state. Software version is displayed: 10 is a software version in the figure.

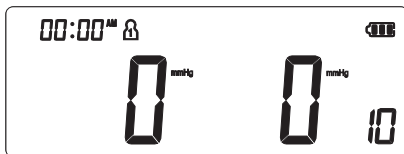**Note:** Only Service personnel is permitted to access this mode. The mode is unavailable in normal use.

## Arm Movement Indicator

If arm movement is detected during the measurement, the “” symbol may appear on the display. This indicates that the measurement result may be inaccurate due to movement. After the measurement is completed, the “” symbol will remain displayed on the LCD screen together with the result.

## Loose Cuff Indicator

When starting the measurement, “” symbol will be displayed if the cuff is applied correctly. If the cuff is wrapped too loosely, the, “” symbol will appear on the display. In this case, reapply the cuff properly and restart the measurement.

## Troubleshooting

| Abnormal phenomenon                              | Cause analysis                                                                                                                        | Processing method                                                                                                                                                              |
|--------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Abnormal operation of the blood pressure monitor | The cuff is tied too tight or too loose, or the cuff is tied incorrectly.                                                             | Adjust the cuff correctly.                                                                                                                                                     |
|                                                  | Moving the arm during measurement or using the electronic sphygmomanometer incorrectly.                                               | Stay quiet, keep your arm steady, and do not move the monitor.                                                                                                                 |
|                                                  | Speaking, nervous or emotional during measurement.                                                                                    | Instead of talking, take deep breaths to calm your mood and relax your body.                                                                                                   |
|                                                  | Incorrect measurement posture.                                                                                                        | Adjust posture, see Page 15.                                                                                                                                                   |
|                                                  | There is interference in charging process or improper operation in measuring process.                                                 | See operation Instructions.                                                                                                                                                    |
|                                                  | If the operator is in a strong electromagnetic interference environment, resulting in the loss or reduction of function/ performance. | It is recommended that the operator, away from the environment, choose an environment with a single electromagnetic field, and then re-measure after the product is recovered. |

The following table shows the error signs that may occur during measurement, possible causes and handling methods. Please measure again using the correct method.


| Error display                                                                    | The cause of the problem                                                                                   | The solution                                                                                                                                                                                      |
|----------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Er1                                                                              | Can't detect high and low Pressure.                                                                        | Please fasten the cuff before measuring.                                                                                                                                                          |
| Er2                                                                              | Cuff too tight or too loose.                                                                               | Please fasten the cuff before measuring.                                                                                                                                                          |
| Er3                                                                              | Arm or body movements or other behaviors may cause abnormal data during the measurement of blood pressure. | Please hold the arm or body still. Please keep quiet and measure again.                                                                                                                           |
| Er4                                                                              | The pressure exceeds 300mmHg.                                                                              | Please fasten the cuff before measuring.                                                                                                                                                          |
| Er5                                                                              | The pressure remains above 15mmHg for 3 minutes.                                                           | Check whether the cuff or air tube is kinked or twisted, and ensure that the air release (deflation) valve is not blocked. If the problem persists, discontinue use and contact the manufacturer. |
| Er6                                                                              | Blood pressure measurements were out of range.                                                             | Please tighten the cuff securely and repeat the measurement. If the problem persists, contact the manufacturer/retailer.                                                                          |
|  | Battery is dead.                                                                                           | Charge the battery by connecting the device to the power adapter (if provided).                                                                                                                   |

**Note:** If you cannot solve the abnormal situation by yourself, you can consult the manufacturer or the manufacturer's designated retailer by phone. It is forbidden to disassemble and repair without permission. If necessary, professional maintenance personnel can ask the manufacturer for the list of components and circuit schematic diagram.

# Blood Pressure Information

## Blood Pressure

Blood pressure is the force of blood pushing against the walls of arteries. It is typically measured in millimeters of mercury (mmHg.) Systolic blood pressure is the maximum force exerted against blood vessel walls each time the heart beats. Diastolic blood pressure is the force exerted on blood vessels when the heart is resting between beats. An individual's blood pressure frequently changes throughout the course of a day. Excitement and tension can cause blood pressure to rise, while drinking alcohol and bathing can lower blood pressure. Certain hormones like adrenaline (which your body releases under stress) can cause blood vessels to constrict, leading to a rise in blood pressure. If these measuring values become too high, it means the heart is working harder than it should.

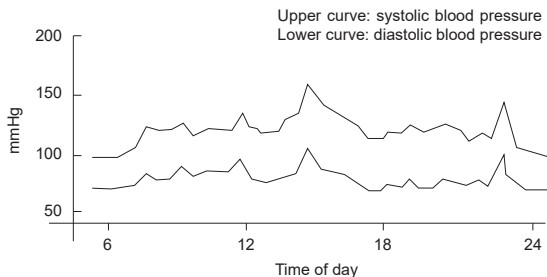
Example: fluctuation within a day (male, 35 years old)

## Operating principles


**Oscillogram method:** Also known as pressure oscillation method, this process is to inflate the cuff to block arterial blood flow, and then detect the gas pressure in the cuff during the process of deflating and extract weak pulse waves. The sphygmomanometer is used to measure blood pressure in the process of pressor, so the control algorithm and structure are needed to make the pressure process stable and uniform. The measurement principle is also oscilloscopic method, which is consistent with the antihypertensive sphygmomanometer. When the measurement results, then open the valve to quickly deflate. The corresponding relationship

between the oscillatory wave amplitude and the pressure waveform in the oscillogram method of the booster electronic sphygmomanometer is shown in the figure below.

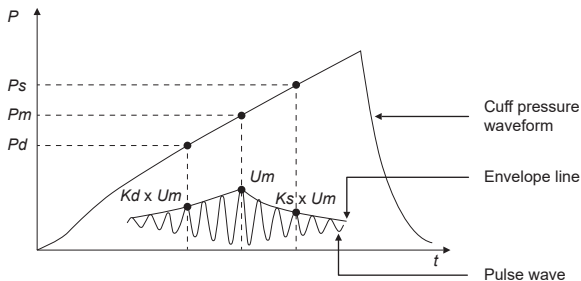
## WHO Blood Pressure Classification Indicator

The Experience M550 AFib is equipped with a classification indicator based on established guidelines from the World Health Organization. The chart below indicates test results (classification indicator is coded also with LCD light on monitor device).

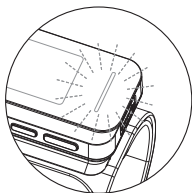
|        |   |                         |
|--------|---|-------------------------|
| Red    | 3 | → Severe Hypertension   |
|        | 2 | → Moderate Hypertension |
| Yellow | 1 | → Mild Hypertension     |
| Green  |   | → High-normal           |
|        | N | → Normal                |
|        |   | → Optimal               |

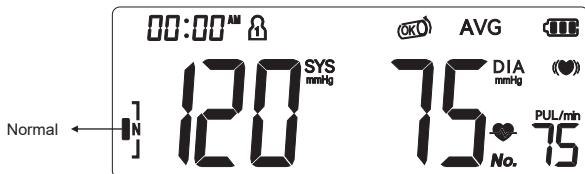**N** : Blood Pressure Classification Indicator

## Health Reminder

Hypertension is a dangerous disease that can affect the quality of life. It can lead to a lot of problems including heart failure, kidney failure, and cerebral hemorrhaging. By maintaining a healthy lifestyle and visiting your physician on a regular basis, hypertension and related diseases are much easier to control when diagnosed in their early stages.

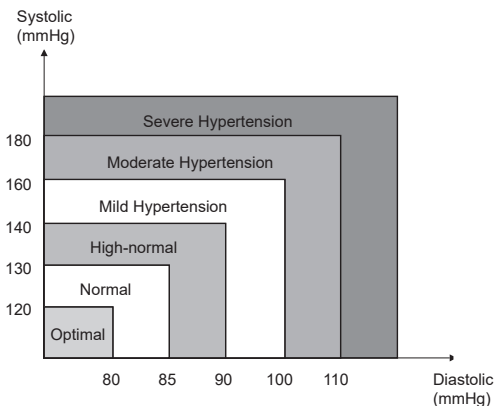**Note:** Do not be alarmed if an abnormal reading occurs. A better indication of an individual's blood pressure occurs after 2–3 readings are taken at the same time each day over an extended period of time. Consult your physician if test results remain abnormal.

## Blood Pressure Q&A

---

**Q:** What is the difference between measuring blood pressure at home or at a professional healthcare clinic?

**A:** Blood pressure readings taken at home are now seen to give a more accurate account as they better reflect your daily life. Readings can be elevated when taken in a clinical or medical environment. This is known as White Coat Hypertension and may be caused by feeling anxious or nervous.

**Note:** Abnormal test results may be caused by:

**1. Improper cuff placement**

Make sure cuff is snug – not too tight or too loose.

Make sure bottom of the cuff is approximately 1–2cm (1/2") above the elbow joint.

**2. Improper body position**

Make sure to keep your body in an upright position.

**3. Feeling anxious or nervous**

Take 2–3 deep breaths, wait a few minutes and resume testing.

**Q:** What causes different readings?

**A:** Blood pressure varies throughout the course of a day. Many factors including diet, stress, cuff placement, etc. may affect an individual's blood pressure.

**Q:** Should I apply the cuff to the left or right arm? What is the difference?

**A:** Either arm can be used when testing, however, when comparing results, the same arm should be used. Testing on your left arm may provide more accurate results as it is located closer to your heart.

**Q:** What is the best time of day for testing?

**A:** Morning time or any time you feel relaxed and stress free.


## Maintenance

---

1. Avoid dropping, slamming, or throwing the device.

2. Avoid extreme temperatures.  
Do not expose the device to direct sunlight.

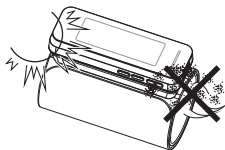3. When cleaning the device, use a soft fabric and lightly wipe with mild detergent. Use a damp cloth to remove dirt and excess detergent.

4. Cuff Cleaning: Do not soak the cuff in water. Apply a small amount of 70% rubbing alcohol to a soft cloth and gently clean the surface of the cuff. Then wipe it with a clean, damp (water-moistened) cloth to remove any residue. Allow the cuff to air dry at room temperature naturally. The cuff must be completely dry before the next use and before being used by another user.

5. Do not use petrol, thinners or similar solvents.

6. Do not disassemble the device.

7. It is recommended that the performance should be checked every 2 years.
8. Expected service life: Approximately three years at 10 tests per day.
9. No service and maintenance while it is in use. Maintenance can be performed only by service personnel. Service and maintenance require parts, repair, technical support will be provided.
10. If it is not used for a long time, it is recommended to charge the lithium battery within 3–4 months to ensure that it does not lose power. Users cannot replace the lithium battery by themselves. If it is confirmed that the lithium battery is faulty, please do not repair it by yourself. Please contact the manufacturer or dealer in time.
11. If the lithium battery is not charged in time, or the lithium battery is placed for a long time, there may be the situation that the lithium battery is about to run out. In this case, the device cannot be turned on. In this case, the lithium battery should be charged immediately.

**The following rules must be observed during use:**

1. Prevent lithium batteries from entering water. Do not heat lithium batteries, do not throw lithium batteries into fire.
2. Lithium batteries should not be used for a long time in an environment of more than 60°C.
3. The device can only use lithium batteries specified by the manufacturer.
4. Use the specified AC adapter to charge the device.
5. To prevent self-discharge during storage, recharge the lithium battery every 3–4 months to maintain normal operation.

## Specifications

|                            |                                                                                                                     |                     |
|----------------------------|---------------------------------------------------------------------------------------------------------------------|---------------------|
| <b>Product Description</b> | Arm-type Fully Automatic Blood Pressure Monitor                                                                     |                     |
| <b>Trade name</b>          | Mediblink Blood pressure monitor Experience M550 AFib                                                               |                     |
| <b>Model/REF</b>           | DBP-61F4L/M550                                                                                                      |                     |
| <b>Display</b>             | LCD Digital Display Size: 79mmx27.5mm                                                                               |                     |
| <b>Measurement Method</b>  | Oscillometric Method                                                                                                |                     |
| <b>Measurement Range</b>   | Systolic Pressure                                                                                                   | 60mmHg~260 mmHg     |
|                            | Diastolic Pressure                                                                                                  | 40mmHg~200 mmHg     |
|                            | Pressure                                                                                                            | 0mmHg~299mmHg       |
|                            | Pressure Accuracy                                                                                                   | ±3mmHg              |
|                            | Pulse                                                                                                               | 30~180 Beats/Minute |
|                            | Pulse Accuracy                                                                                                      | ±5%                 |
| <b>Pressurization</b>      | Automatic Pressurization                                                                                            |                     |
| <b>Memory</b>              | 2x150 Memories in Two Groups with Date and Time                                                                     |                     |
| <b>Function</b>            | Irregular Heartbeat Detection                                                                                       |                     |
|                            | WHO Classification Indicator (graph and LED indicator)                                                              |                     |
|                            | Last 3 Tests Average                                                                                                |                     |
|                            | Low Battery Detection                                                                                               |                     |
|                            | Automatic Power-Off                                                                                                 |                     |
|                            | Backlight                                                                                                           |                     |
|                            | Mean value measurement function                                                                                     |                     |
| <b>Power Source</b>        | AFib detection function                                                                                             |                     |
|                            | Lithium battery 3.7 V or DC 5.0V,1000mA medical external power supply (DC 5.0V, 1000mA) (recommended, not provided) |                     |

|                                    |                                                                                                                                                   |                          |
|------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------|
| <b>Battery Life</b>                | Approximately 2 months at 3 tests per day                                                                                                         |                          |
| <b>Device Weight</b>               | Approximately 241 g (8.5 oz) (excluding battery)                                                                                                  |                          |
| <b>Device Dimensions</b>           | Approx. 130mmx49mmx25.5mm (5.12"x1.93"x1.00") (LxWxH)                                                                                             |                          |
| <b>Cuff Circumference</b>          | Pre-formed/Rigid arm cuff (Medium):<br>Fits arm circumference 22–42cm.                                                                            |                          |
| <b>Operating Environment</b>       | Temperature                                                                                                                                       | 10°C~40°C (50°F~104°F)   |
|                                    | Humidity                                                                                                                                          | 15%~93%RH                |
|                                    | Pressure                                                                                                                                          | 80kPa~106kPa             |
| <b>Storage Environment</b>         | Temperature                                                                                                                                       | -25°C~55°C (-13°F~131°F) |
|                                    | Humidity                                                                                                                                          | ≤93%RH                   |
| <b>Classification</b>              | Internal Powered Equipment, Type BF <br>Cuff is the Applied Part |                          |
| <b>Ingress Protection rating</b>   | IP 22, Indoor Used Only                                                                                                                           |                          |
| <b>Battery Storage Temperature</b> | -25°C~55°C (-13°F~131°F)                                                                                                                          |                          |

This Blood Pressure Monitor complies with the European regulations and bears the CE mark “CE 0123”. This blood pressure monitor also complies with mainly following standards (included but not limited):

**Safety standard:**

EN 60601-1 Medical electrical equipment part 1: General requirements for safety EMC standard:

EN 60601-1-2 Medical Electrical Equipment – Part 1–2: General Requirements for Basic Safety and Essential Performance – Collateral Standard: Electromagnetic Disturbances – Requirements and Tests.

**Performance standards:**

IEC80601-2-30, Medical electrical equipment – Part 2–30: Particular requirements for the basic safety and essential performance of automated non-invasive sphygmomanometers.

ISO 81060-2, non-invasive sphygmomanometers – part 2: clinical validation of automated measurement type.

**Other standards:**

EN ISO 15223-1, Medical devices – Symbols to be used with information to be supplied by the manufacturer – Part 1: General requirements

## Electromagnetic Compatibility Information

The device complies with EMC requirements of the international standard IEC 60601-1-2. Compliance is demonstrated under the conditions specified in the tables below. The device is an electrical medical product and is subject to special precautionary measures with regard to EMC which must be published in the instructions for use. Portable and mobile HF communications equipment can affect the device. Use of the device in conjunction with non-approved accessories can affect the device negatively and alter the electromagnetic compatibility. Environment of Intended use: Home healthcare environment and Professional healthcare facility environment.

**Table 1:**

| Guidance and declaration of manufacturer-electromagnetic emissions                                                                                                              |                 |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------|
| The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment. |                 |
| Emissions test                                                                                                                                                                  | Compliance      |
| Radiated emission CISPR 11                                                                                                                                                      | Group 1,class B |
| Conducted emission CISPR 11                                                                                                                                                     | Group 1,class B |
| Harmonic emissions IEC 61000-3-2                                                                                                                                                | Class A         |
| Voltage fluctuations/flicker emissions IEC 61000-3-3                                                                                                                            | Complies        |

Table 2:

| The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment. |                                                                                                                                                                 |                                                                                                                                                                |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------|
| IMMUNITY test                                                                                                                                                                   | IEC 60601 test level                                                                                                                                            | Compliance level                                                                                                                                               |
| Electrostatic discharge (ESD)<br>IEC 61000-4-2                                                                                                                                  | $\pm 8$ kV<br>contact discharge<br>$\pm 2$ kV, $\pm 4$ kV, $\pm 8$ kV, $\pm 15$ kV air                                                                          | $\pm 8$ kV<br>contact discharge<br>$\pm 2$ kV, $\pm 4$ kV, $\pm 8$ kV, $\pm 15$ kV air                                                                         |
| Electrical fast transients / bursts<br>IEC 61000-4-4                                                                                                                            | $\pm 2$ kV, 100 kHz, for AC power port                                                                                                                          | $\pm 2$ kV, 100 kHz, for AC power port                                                                                                                         |
| Surges<br>IEC 61000-4-5                                                                                                                                                         | $\pm 0.5$ kV, $\pm 1$ kV (differential mode)                                                                                                                    | $\pm 0.5$ kV, $\pm 1$ kV (differential mode)                                                                                                                   |
| Voltage dips, short interruptions<br>IEC 61000-4-11                                                                                                                             | 0 % UT; 0.5 cycle<br>At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% UT; 1 cycle and<br>70% UT; 25/30 cycles<br>Single phase: at 0°<br>0% UT; 250/300 cycle | 0% UT; 0.5 cycle<br>At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% UT; 1 cycle and<br>70% UT; 25/30 cycles<br>Single phase: at 0°<br>0% UT; 250/300 cycle |
| Rated power frequency (50/60Hz)<br>magnetic fields<br>IEC 61000-4-8                                                                                                             | 30 A/m; 50Hz or 60Hz                                                                                                                                            | 30 A/m; 50Hz or 60Hz                                                                                                                                           |
| Radiated RF EM fields<br>IEC 61000-4-3                                                                                                                                          | 3V/m and 10 V/m<br>80MHz–2.7GHz<br>80%AM at 1kHz                                                                                                                | 3V/m and 10 V/m<br>80MHz–2.7GHz<br>80%AM at 1kHz                                                                                                               |
| Proximity fields from RF wireless communications equipment<br>IEC 61000-4-3                                                                                                     | See Table 3                                                                                                                                                     | See Table3                                                                                                                                                     |
| Conducted disturbances induced by RF fields<br>IEC 61000-4-6                                                                                                                    | 3 V in 0.15MHz – 80MHz<br>6 V in ISM and amateur radio bands between 0.15MHz and 80MHz<br>80% AM at 1kHz                                                        | 3 V in 0.15MHz – 80MHz<br>6 V in ISM and amateur radio bands between 0.15MHz and 80MHz<br>80%AM at 1kHz                                                        |
| Proximity magnetic fields<br>IEC 61000-4-39                                                                                                                                     | See Table 4                                                                                                                                                     | See Table 4                                                                                                                                                    |

Table 3:

| Guidance and declaration of manufacturer-electromagnetic immunity                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |            |                                                                          |                                     |                           |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------|--------------------------------------------------------------------------|-------------------------------------|---------------------------|
| Nowadays, many RF wireless equipment have been used in various healthcare locations where medical equipment and/or systems are used. When they are used in close proximity to medical equipment and/or systems, the medical equipment and/or systems' basic safety and essential performance may be affected. Arm-type Fully Automatic Digital Blood Pressure Monitor has been tested with the immunity test level in the below table and meet the related requirements of IEC 60601-1-2:2014. The customer and/or user should help keep a minimum distance between RF wireless communications equipment and this medical equipment and/or systems as recommended below. |            |                                                                          |                                     |                           |
| Test frequency (MHz)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | Band (MHz) | Service                                                                  | Modulation                          | Immunity test level (V/m) |
| 385                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | 380–390    | TETRA 400                                                                | Pulse modulation 18Hz               | 27                        |
| 450                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | 430–470    | GMRS 460<br>FRS 460                                                      | FM<br>± 5kHz deviation<br>1kHz sine | 28                        |
| 710                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | 704–787    | LTE Band 13, 17                                                          | Pulse modulation 217Hz              | 9                         |
| 745                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |            |                                                                          |                                     |                           |
| 780                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |            |                                                                          |                                     |                           |
| 810                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | 800–960    | GSM 800/900,<br>TETRA 800, iDEN<br>820, CDMA 850, LTE<br>Band 5          | Pulse modulation 18Hz               | 28                        |
| 870                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |            |                                                                          |                                     |                           |
| 930                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |            |                                                                          |                                     |                           |
| 1720                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 1700–1990  | GSM 1800; CDMA<br>1900; GSM 1900;<br>DECT; LTE Band 1,<br>3, 4, 25; UMTS | Pulse modulation 217Hz              | 28                        |
| 1845                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |            |                                                                          |                                     |                           |
| 1970                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |            |                                                                          |                                     |                           |
| 2450                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 2400–2570  | WLAN, 802.11 b/g/n,<br>RFID 2450, LTE<br>Band 7                          | Pulse modulation 217Hz              | 28                        |
| 5240                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 5100–5800  | WLAN 802.11<br>a/n                                                       | Pulse modulation 217Hz              | 9                         |
| 5500                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |            |                                                                          |                                     |                           |
| 5785                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |            |                                                                          |                                     |                           |

Table 4:

| Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields |                                          |                          |
|------------------------------------------------------------------------------|------------------------------------------|--------------------------|
| Test frequency                                                               | Modulation                               | IMMUNITY TEST LEVEL(A/m) |
| 30kHz <sup>a)</sup>                                                          | CW                                       | 8                        |
| 134,2kHz                                                                     | Pulse Modulation <sup>b)</sup><br>2,1KHz | 65 <sup>c)</sup>         |
| 13,56MHz                                                                     | Pulse Modulation <sup>b)</sup><br>50KHz  | 7,5 <sup>d)</sup>        |

<sup>a)</sup> This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME HEALTHCARE ENVIRONMENT.

<sup>b)</sup> The carrier shall be modulated using a 50% duty cycle square wave signal.

<sup>c)</sup> r.m.s., before modulation is applied.

## Additional Notes

---

1. **WARNING:** Use of this equipment adjacent to or stacked with other equipment should be avoided, as it may result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
2. **WARNING:** Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Arm-Type Fully Automatic Digital Blood Pressure Monitor, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment may occur.
3. The software identifier refers to the software evaluation report.
4. **Verification of manometer pressure accuracy:** In the power-off state, press and hold the "START/STOP" button, then install the batteries. When the LCD screen is fully displayed, release the "START/STOP" button. When the LCD shows "00", the device enters static mode. At this point, a 500 ml air capacity device, a calibrated standard pressure gauge, and a manual pressure device can be connected to the sphygmomanometer via the cuff interface. Apply pressure manually within the effective measurement range and compare the readings of the sphygmomanometer with those of the standard pressure gauge. This mode can be used to verify manometer pressure accuracy.
5. **Contraindications:** This product is not intended for use by infants or individuals who are unable to express their intentions.
6. **Intended Use:** This device is intended for home use and professional healthcare use.
7. **Patient as operator:** The patient is an intended operator. The patient must not perform maintenance operations other than battery replacement or recharging (if applicable).
8. **WARNING:** Do not modify this equipment without authorization from the manufacturer.
9. **Essential performance – maintenance advice:** Pressure calibration is performed before the product leaves the factory. Users may verify accuracy using the method described in the section "Verification of Manometer Pressure Accuracy". If significant deviation is observed, contact the manufacturer for recalibration.
10. **Mechanical strength and heat resistance:** The device will maintain its mechanical strength and resistance to heat throughout its expected service life.
11. The cuff and housing of the blood pressure monitor have been tested for biocompatibility and do not contain allergenic or harmful materials. Discontinue use if an allergic reaction occurs.

12. **WARNING:** Non-professionals must not modify the equipment, as this may result in inaccurate measurements.
13. **WARNING:** Do not expose the device to prolonged environmental stress (e.g. direct sunlight, high temperature, or humidity), as this may reduce its performance.
14. **WARNING:** This device must not be used by children or animals.
15. **Cleaning:** The device can be cleaned by the lay operator according to the cleaning procedures described in this manual.
16. **WARNING:** Do not use a damaged cuff for blood pressure measurement.
17. **WARNING:** If the user experiences severe discomfort during measurement, press the button to deflate the cuff immediately or remove the cuff from the arm.
18. **WARNING:** If an unexpected reading occurs, the operator may take additional measurements and consult a healthcare professional.
19. **WARNING:** Use of the device outside the specified environmental conditions may damage the device and result in inaccurate measurements.
20. This ME equipment is not intended for use in the presence of flammable anesthetic mixtures or in oxygen-rich environments.
21. The operator shall ensure that basic safety and essential performance are maintained with regard to electromagnetic compatibility (EMC).
22. **WARNING:** The operator must discontinue use and contact customer service if essential performance is lost or degraded due to electromagnetic disturbances.
23. **WARNING:** Use of accessories, transducers, or cables other than those specified or provided by the manufacturer may result in increased electromagnetic emissions or decreased electromagnetic immunity, leading to improper operation.

### **Correct Disposal of This Product** (Waste Electrical & Electronic Equipment)

This marking shown on the product indicates that it should not be disposed with other household waste at the end of its life. To prevent potential harm to the environment or to human health, please separate this product from other types of wastes and recycle it responsibly. When disposing of this type of product, contact the retailer where product was purchased or contact your local government office for details regarding how this item can be disposed in an environmentally safe recycling center.

Business users should contact their supplier and check the terms and conditions of the purchasing agreement. This product should not be mixed with other commercial wastes for disposal. This product is free of hazardous materials.

# Warranty

**Product: Mediblink Blood Pressure Monitor Experience M550 AFib**

**Manufactured for (Importer for EU & Distributor)** Mediblink d.o.o.,  
Gubčeva cesta 19, 8210 Trebnje, Slovenia; [info@mediblink.com](mailto:info@mediblink.com); [www.mediblink.com](http://www.mediblink.com)

**Seller's name, address, signature and stamp\*:**

**Date of delivery/sales\*:**

\*If the invoice is accompanied by this warranty, and if all above information can be seen from the invoice, it is not necessary to fill in this field.

## WARRANTY TERMS

Dear customers!

The warranty period for the main unit is **10 years** and starts from the date of purchase or delivery of the product to the customer. The warranty period for the cuff and the rechargeable battery is **2 years** and starts from the date of purchase or delivery of the product to the customer. The cuff and the rechargeable battery are considered wearable components with a limited service life. The warranty applies to original products purchased from sellers or distributors who provide valid proof of purchase and product traceability. When claiming the warranty, the customer must provide the invoice or other proof of purchase, or a completed warranty certificate. The consumer's statutory rights arising from the conformity of goods remain unaffected.

Many complaints are the result of improper use of the product. You can easily avoid any problem, by getting useful information provided by our special service department. To reach our service department, you can call or send e-mail to Mediblink local distributor.

Before sending the product back to retailer, we kindly ask you to call our service department, to get help about how to use the device to save you with unneeded trips.

The manufacturer guarantees free elimination of all imperfections due to defects in material or manufacturing procedure by repairing or replacing the product. In case that the product can not be repaired or replaced, the customer will get the money refund. The guarantee is not valid in case of the force majeure, accidents or unexpected events (such as lightning, water, fire etc.), incorrect use or incorrect transport, non-compliance with safety and maintaining regulations or in case of unprofessional product intervention.

Traces of every day product usage (scratches, abrasions) and not subject to claim. The warranty does not eliminate the customer rights, which originate from seller responsibility for product flaws. By accepting the claimed product by the service department, the service department does not take responsibility for loss of saved data or settings on the product. All product repairs, which are performed out of product warranty period, have to be paid by customer by prior notice.

The manufacturer guarantees the product quality and flawless product operation in the warranty period, which starts with the day of product purchase. If the product can not be repaired in 45 days, the product will be replaced with a new one. In case that the product can not be replaced, the money will be refunded to the customer.

In case of product claim, call or send e-mail to Mediblink local distributor.

**Replacement parts can be ordered at [www.mediblink.com](http://www.mediblink.com).**

**Manufacturer:**  
JOYTECH Healthcare Co., Ltd.  
No.365, Wuzhou Road, 311100  
Hangzhou, Zhejiang Province,  
People's Republic of China  
E-mail: [info@sejoy.com](mailto:info@sejoy.com)

**Manufactured for**  
(Importer for EU & Distributor):  
Mediblink d.o.o.,  
Gubčeva cesta 19,  
8210 Trebnje,  
Slovenia  
[info@mediblink.com](mailto:info@mediblink.com);  
[www.mediblink.com](http://www.mediblink.com)

**Shanghai International Holding**  
Corp. GmbH (Europe)  
Eiffestrasse 80, 20537 Hamburg,  
Germany  
E-mail: [shholding@hotmail.com](mailto:shholding@hotmail.com)

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