**EN Mediblink Wrist Blood Pressure Monitor M505**

**INSTRUCTIONS FOR USE**

Please read the instructions for use carefully before using the product

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Pažljivo pročitajte upute prije upotrebe uređaja

# TABLE OF CONTENTS

|                                                                                      |           |
|--------------------------------------------------------------------------------------|-----------|
| <b>1. Introduction</b>                                                               | <b>6</b>  |
| 1.1 Features of M505                                                                 | 6         |
| 1.2 Important information about self-measurement                                     | 6         |
| <b>2. Important Information on the Subject of Blood-Pressure and Its Measurement</b> | <b>7</b>  |
| 2.1 How does high/low blood-pressure arise?                                          | 7         |
| 2.2 Which values are normal?                                                         | 7         |
| <b>3. The Various Components of the Blood-Pressure Monitor</b>                       | <b>9</b>  |
| <b>4. Putting the Blood-Pressure Monitor Into Operation</b>                          | <b>10</b> |
| 4.1 Inserting the batteries                                                          | 10        |
| 4.2 Reading the set date                                                             | 10        |
| <b>5. Carrying Out a Measurement</b>                                                 | <b>14</b> |
| 5.1 Before the measurement                                                           | 14        |
| 5.2 Common sources of error                                                          | 14        |
| 5.3 Fitting the cuff                                                                 | 15        |
| 5.4 Measuring procedure                                                              | 15        |
| 5.5 Discontinuing a measurement                                                      | 16        |
| 5.6 Memory – storage and recall of the measurements                                  | 16        |
| 5.7 Memory – cancellation of all measurements                                        | 16        |
| <b>6. Error Messages/Malfunctions</b>                                                | <b>17</b> |
| <b>7. Care and Maintenance, Recalibration</b>                                        | <b>19</b> |
| <b>8. Service Life</b>                                                               | <b>19</b> |
| <b>9. Battery Life</b>                                                               | <b>19</b> |
| <b>10. Safety, Care and Disposal</b>                                                 | <b>20</b> |
| <b>11. Reference to Standards</b>                                                    | <b>21</b> |
| <b>12. Technical Specifications</b>                                                  | <b>21</b> |
| <b>13. Manufacturer's Declaration</b>                                                | <b>23</b> |
| <b>14. Warranty</b>                                                                  | <b>27</b> |

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---

## 1.1 Features of M505

The blood-pressure monitor Wrist Blood Pressure Monitor M505 with MDI, measuring during inflation (with integrated time/date display) is a fully automatic, digital blood-pressure measuring device for use on the wrist, which enables very fast and reliable measurement of the systolic and diastolic blood-pressure as well as the pulse frequency by way of the oscillometric method of measuring. The device offers a very high and clinical tested measurement accuracy and has been designed to provide a maximum of user-friendliness. The device is intended for self-use at home.

Before using, please read through this instruction manual carefully and then keep it in a safe place. For further questions on the subject of blood-pressure and its measurement, please contact your doctor.

### **Warning!**

## 1.2 Important information about self-measurement

- Substitution of a different component might result in measurement error.
- Do not use with neonatal patients.
- Too frequent measurements can cause injury to the patient due to blood flow interference.
- The application of the cuff over a wound can cause further injury.
- The application of the cuff and its pressurization on any limb where intravascular access or therapy, or an arteriovenous (A-V) shunt, is present because of temporary interference to blood flow and could result in injury to the patient.
- Do not let the cuff and its pressurization on the arm on the side of a mastectomy
- The need to check that operation of the automated sphygmomanometer does not result in prolonged impairment of patient blood circulation.
- Not intended to be used together with HF surgical equipment.
- Do not forget: self-measurement means control, not diagnosis or treatment. Unusual values must always be discussed with your doctor. Under no circumstances should you alter the dosages of any drugs prescribed by your doctor.
- The pulse display is not suitable for checking the frequency of heart pacemakers!

- In cases of cardiac irregularity (Arrhythmia), measurements made with this instrument should only be evaluated after consultation with the doctor.

### **Electromagnetic interference**

The device contains sensitive electronic components (Microcomputer). Therefore, avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g. mobile telephones, microwave cookers). These can lead to temporary impairment of the measuring accuracy.

## **2. Important Information on the Subject of Blood-Pressure and Its Measurement**

---

### **2.1 How does high/low blood-pressure arise?**

The level of blood-pressure is determined in a part of the brain, the so-called circulatory center, and adapted to the respective situation by way of feedback via the nervous system. To adjust the blood-pressure, the strength and frequency of the heart (Pulse), as well as the width of circulatory blood vessels is altered. The latter is effected by way of fine muscles in the blood-vessel walls. The level of arterial blood-pressure changes periodically during the heart activity: During the »blood ejection« (Systole) the value is maximal (systolic blood-pressure value), at the end of the heart's »rest period« (Diastole) minimal (diastolic blood-pressure value). The blood-pressure values must lie within certain normal ranges in order to prevent particular diseases.

### **2.2 Which values are normal?**

Blood pressure is too high if at rest, the diastolic pressure is above 90 mmHg and/or the systolic blood-pressure is over 160 mmHg. In this case, please consult your doctor immediately. Long-term values at this level endanger your health due to the associated advancing damage to the blood vessels in your body.

Should the systolic blood-pressure values lie between 140 mmHg and 160 mmHg and/or the diastolic blood-pressure values lie between 90 mmHg and 100 mmHg, likewise, please consult your doctor.

Furthermore, regular self-checks will be necessary.

With blood-pressure values that are too low, i.e. systolic values under 100 mmHg and/or diastolic values under 60 mmHg, likewise, please consult your doctor.

Even with normal blood-pressure values, a regular self-check with your blood-pressure monitor is recommended. In this way you can detect possible changes in your values early and react appropriately. If you are undergoing medical treatment to control your blood pressure, please keep a record of the level of your blood pressure by carrying out regular self-measurements at specific times of the day. Show these values to your doctor.

**Never use the results of your measurements to alter independently the drug doses prescribed by your doctor.**

Table for classifying blood-pressure values (unit: mmHg) according to World Health Organization:


| Range                           | Systolic Blood-pressure | Diastolic Blood-pressure | Measures                      |
|---------------------------------|-------------------------|--------------------------|-------------------------------|
| Blood pressure optimum          | between 100 and 120     | between 60 and 80        | Self-check                    |
| Blood pressure normal           | between 120 and 129     | between 80 and 84        | Self-check                    |
| Blood pressure slightly high    | between 130 and 139     | between 85 and 89        | Consult your doctor           |
| Blood pressure too high         | between 140 and 159     | Between 90 and 99        | Seek medical advice           |
| Blood pressure far too high     | between 160 and 179     | Between 100 and 109      | Seek medical advice           |
| Blood pressure dangerously high | Higher than 180         | Higher than 110          | Urgently seek medical advice! |

### Further information

- If your values are mostly standard under resting conditions but exceptionally high under conditions of physical or psychological stress, it is possible that you are suffering from so-called »labile hypertension«. Please consult your doctor if you suspect that this might be the case.
- Correctly measured diastolic blood-pressure values above 120 mmHg require immediate medical treatment.

### 3. The Various Components of the Blood-Pressure Monitor
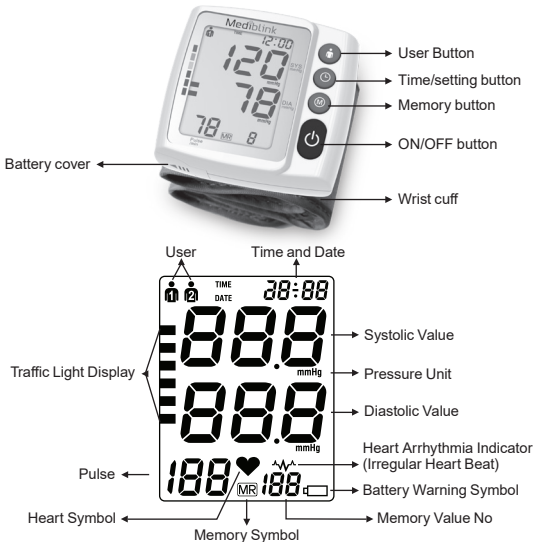
## 4. Putting the Blood-Pressure Monitor Into Operation

### 4.1 Inserting the batteries

- Insert the batteries (2 x size AAA 1.5V), thereby observing the indicated polarity.
- If the battery warning »« appears in the display, the batteries are empty and must be replaced by new ones.

#### **Warning!**

- After the battery warning »« appears, the device is blocked until the batteries have been replaced.
- Please use »AAA« Long-Life or Alkaline 1.5V Batteries. The use of 1.2V Accumulators is not recommended.
- If the blood-pressure monitor is left unused for long periods, please remove the batteries from the device.


### 4.2. Reading the set date

Please press the TIME button, the date will be shown in the display.


#### **User selection:**

This advanced blood pressure monitor allows you to track blood pressure readings for 2 individuals independently.

- Before measurement, make sure you set the unit for the intended user. The unit can track results for 2 individuals (User 1, User 2).
- Press the TIME button for at least 3 seconds. The display now indicates the set user, during which the set user blink. To confirm, press ON/OFF button.

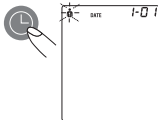- c) Click the MEMORY button to select User.
- d) We suggest the first person to take their pressure to be User 1.

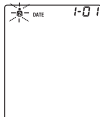
### Setting the time, date

This blood-pressure monitor incorporates an integrated clock with date display. This has the advantage, that at each measurement procedure, not only the blood-pressure values are stored, but also the exact moment of the measurement.


After new batteries have been inserted, the clock begins to run from the following setting: 2018-01-01 12:00 O'clock.

You must then re-enter the date and current time. For this, please proceed as follows:


- 1) Press the TIME button for at least 3 seconds firstly, user icon blinks. Then press TIME button again, and the display now indicates the set year, during which the four characters blink.

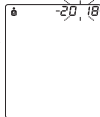
- 2) The correct year can be entered by pressing the MEMORY button.

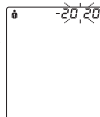
- 3) Press the TIME button again. The display now switches to the current date, during which the first character (month) blinks.

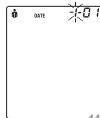
4) The corresponding month can now be entered by pressing the MEMORY button.

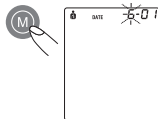
5) Press the TIME button again. The last two characters (day) are now blinking.

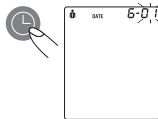
6) The corresponding day can now be entered by pressing the MEMORY button.

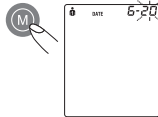
7) Press the TIME button again. The display now switches to the current time, during which the first character (Hour) blinks.

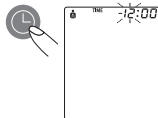
8) The corresponding hour can now be entered by pressing the MEMORY button.

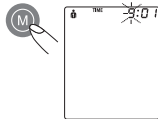
9) Press the TIME button again. The last two characters (Minutes) now blink.

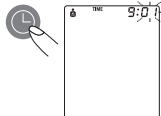10) The exact time can now be entered by pressing the MEMORY button.

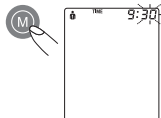11) Now after all settings have been made, press the TIME button once again. The date is briefly displayed and then the time. The input is now confirmed and the clock begins to run.
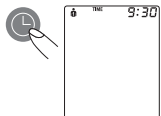
### Further Information

With each press of the button (TIME, MEMORY) one input is made (e.g. switching over from hours to minutes mode, or altering the value by +1).

However, if you keep the respective button depressed, you can switch more quickly to find the desired value respectively.

## 5. Carrying Out a Measurement

---

### 5.1 Before the measurement:

- Avoid eating, smoking as well as all forms of exertion directly before the measurement. All these factors influence the measurement result. Try and find time to relax by sitting in an armchair in a quite atmosphere for about ten minutes before the measurement.
- Measure always on the same wrist (normally left).
- Attempt to carry out the measurements regularly at the same time of day, since the blood-pressure changes during the course of the day.


### 5.2 Common sources of error:

**Comparable blood-pressure measurements always require the same conditions! These are normally always quiet conditions.**

- All efforts by the patient to support the arm can increase the blood-pressure. Make sure you are in a comfortable, relaxed position and do not activate any of the muscles in the measurement. Use a cushion for support if necessary.
- The performance of the automated sphygmomanometer can be affected by extremes of temperature, humidity and altitude.
- If the wrist artery lies considerably lower (higher) than the heart, an erroneously higher (lower) blood-pressure will be measured! (Each 15cm difference in height results in a measurement error of 10 mmHg!)
- A loose cuff causes false measurement values.
- With repeated measurements, blood accumulates in the respective wrist which can lead to false results. Correctly executed blood-pressure measurements should therefore first be repeated after a 5 minute pause or after the arm has been held up in order to allow the accumulated blood to flow away (after at least 3 minutes).

## 5.3 Fitting the cuff

- Remove all eventual objects and jewellery (e.g. wristwatch) from the wrist in question. Draw the cuff over the wrist.
- The distance between the cuff and the hand should be approx. 10 mm
- Secure the cuff with the Velcro fastener, so that it lies comfortably and not too tight, whereby no space should remain between the cuff and the wrist.
- Sit up straight, hold the left arm in front of the chest, with the palm inward. The cuff should be at the same height as the heart. Take care that the cuff lies free. Remain so for 2 minutes sitting quietly, before beginning with the measurement.
- Let legs uncrossed, feet flat on the floor, back and arm supported.
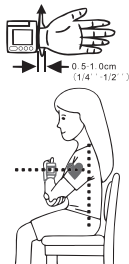
## 5.4 Measuring procedure

After the cuff has been appropriately positioned, the measurement can begin:

- Press the ON/OFF button, the pump begins to inflate the cuff. In the display, the increasing cuff-pressure is continually displayed.
- As the cuff inflates, the monitor automatically determines your ideal inflation level. This monitor detects your blood pressure and pulse rate during inflation. The heartbeat symbol (♥) flashes at every heartbeat.
- When the measurement has been concluded, The measured systolic and diastolic blood-pressure values as well as the pulse frequency are now displayed.

Example (Fig.): Systole 118, Diastole 73, Pulse 75

The measurement results are displayed, until you switch the device off. If no button is pressed for 3 minutes, the device switches automatically off, to save the batteries.

## 5.5 Discontinuing a measurement

If it is necessary to interrupt a blood pressure measurement for any reason (e.g. the patient feels unwell), the »ON/OFF« button can be pressed at any time. The device then immediately lowers the cuff-pressure automatically.

## 5.6 Memory – storage and recall of the measurements

The blood-pressure monitor automatically stores each of the last 120 measurement values. By pressing the MEMORY button, **an average value** of the last 3 measurements as well as the last measurement (MR1) and the further last 119 measurements (MR2, MR3 ..., MR120) can be displayed one after the other.


(MR1: Values of the last measurement)

(MR2-MR120: Values of the measurement before MR1)

## 5.7 Memory– cancellation of all measurements Attention!

Before you delete all readings stored in the memory, make sure you will not need refer to the readings at a later date. Keeping a written record is prudent and may provide additional information for your doctor's visit. In order to delete all stored readings, depress the MEMORY button for at least 5 seconds, the display will show the symbol »CL« and then release the button. to permanently clear the memory, Press the MEMORY button while »CL« is flashing.
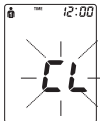
## 6. Error Messages/Malfunctions

If an error occurs during a measurement, the measurement is discontinued and a corresponding error code is displayed (Example: Error No. 2).

| Error No. | Possible cause(s)                                                                                                                                                                                                      |
|-----------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| ERR 1     | No pulse has been detected.                                                                                                                                                                                            |
| ERR 2     | Unnatural pressure impulses influence the measurement result. Reason: hand was moved during the Measurement (Artefact).                                                                                                |
| ERR 3     | 1) The inflation takes too long. The cuff is not correctly seated.<br>2) Positioning Detection error.                                                                                                                  |
| ERR 5     | The measured readings indicated an unacceptable difference between systolic and diastolic pressures. Take another reading following directions carefully. Contact your doctor if you continue to get unusual readings. |
| ERR 8     | Pressure in cuff is over 290 mmHg.                                                                                                                                                                                     |

### Other possible malfunctions and their elimination

If problems occur when using the device, the following points should be checked and if necessary, the corresponding measures are to be taken:

| Malfunction                                                                                       | Remedy                                                                                                                                                                                                    |
|---------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| The display remains empty when the instrument is switched on although the batteries are in place. | <ol style="list-style-type: none"> <li>1. Check batteries for correct polarity and if necessary insert correctly.</li> <li>2. If the display is unusual, re-insert batteries or exchange them.</li> </ol> |

|                                                                                                                                     |                                                                                                                                                                                                                                                                                                                     |
|-------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <p>The device frequently fails to measure the blood pressure values, or the values measured are too low (too high).</p>             | <ol style="list-style-type: none"> <li>1. Check the positioning of the cuff.</li> <li>2. Measure the blood-pressure again in peace and quiet under observance of the details made under point 5.</li> </ol>                                                                                                         |
| <p>Every measurement produces a different value although the instrument functions normally and the values displayed are normal.</p> | <ol style="list-style-type: none"> <li>1. Please read the following information and the points listed under »Common sources of error«. Repeat the measurement.<br/><b>Please note: Blood pressure fluctuates continually so successive measurements will show some variability.</b></li> </ol>                      |
| <p>Blood pressure measured differs from those values measured by the doctor.</p>                                                    | <ol style="list-style-type: none"> <li>1. Record the daily development of the values and consult your doctor.<br/><b>Please note: Individuals visiting their doctor frequently experience anxiety which can result in a higher reading at the doctor than obtained at home under resting conditions.</b></li> </ol> |

## Further Information

The level of blood-pressure is subject to fluctuations even with healthy people. Important thereby is, that comparable **measurements always require the same conditions (Quiet conditions)!**

If, in spite of observing all these factors, the fluctuations are larger than 15 mmHg, and/or the arrhythmia symbol  is in the display, please consult your doctor. For licensing, the device has been subjected to strict clinical tests, by which the computer program used to measure the blood-pressure values was tested by experienced specialist doctors in Germany. The same computer program is used in every individual device, and has thus also been clinically tested. The manufacture of the devices takes place according to the terms of the European standard for blood-pressure measuring devices. You must consult your local distributor if there are technical problems with the blood-pressure instrument. **Never attempt to repair the instrument yourself!**

Any unauthorised opening of the instrument invalidates all guarantee claims!

## 7. Care and Maintenance, Recalibration

---

- a) Do not expose the device to either extreme temperatures, humidity, dust or direct sunlight.
- b) The cuff contains a sensitive air-tight bubble. Handle this carefully and avoid all types of straining through twisting or buckling.
- c) Clean the device with a soft, dry cloth. Do not use petrol, thinners or similar solvent. Spots on the cuff can be removed carefully with a damp cloth and soapsuds. The cuff must not be washed!
- d) Do not drop the instrument or treat it roughly in any way. Avoid strong vibrations.
- e) Never open the device! Otherwise the manufacturer calibration becomes invalid!

## 8. Service Life

---

5 years.

## 9. Battery Life

---

300 times measurement with 2- size »AAA« alkaline Batteries.

## 10. Safety, Care and Disposal

### **Safety and protection**

- This instrument may be used only for the purpose described in this booklet. The manufacturer cannot be held liable for the damage caused by incorrect application.
- This instrument comprise sensitive components and must be treated with caution. Observe the storage and operating condition described in the »Technical specifications« section!
- Protect it from:
  - water and moisture
  - extreme temperatures
  - impact and dropping
  - contamination and dust
  - direct sunlight
  - heat and cold
- The cuffs are sensitive and must be handled with care.
- Only pump up the cuff once fitted.
- Do not use the instrument close to strong electromagnetic fields such as mobile telephones or radio installations.
- Do not use the instrument if you think it is damaged or notice anything unusual.
- If the instrument is not going to be used for a prolonged period the batteries should be removed.
- Read the additional safety instructions in the individual sections of this booklet.
-  Ensure that children do not use the instrument unsupervised: some parts are small enough to be swallowed.
- Must use the recognized accessories, detachable parts and materials, if the use of other parts or materials can degrade minimum safety.
- A warning to remove primary batteries if the instruments is not likely to be used for some time.

### **Instrument care**

Clean the instrument only with a soft, dry cloth.

### **Disposal**

Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

## 11. Reference to Standards

Device standards: Device corresponds to the requirements of the European standard for non-invasive blood- pressure monitor.

### Standards

EN IEC 80601-2-30, EN ISO 81060-2, EN 60601-1, EN 60601-1-6, EN 62366-1, EN 62304, EN ISO 10993-1, EN ISO 10993-5, EN ISO 10993-23, ISO 14971, EN 60601-1-2, EN 60601-1-11

The Device meets the requirements of Directive 93/42/EEC for Medical Devices Class IIa.

## 12. Technical Specifications

|                        |                                                                                            |
|------------------------|--------------------------------------------------------------------------------------------|
| Model / REF:           | BP606W / M505                                                                              |
| Measurement Procedure: | Oscillometric, corresponding to Korotkoff method:<br>Phase I: systolic, Phase V: diastolic |
| Display:               | Digital display                                                                            |
| Measuring range:       | Pressure: 30 to 280 mmHg (in 1 mmHg increment)<br>Pulse: 40 to 199 beat/minute             |
| Static accuracy:       | Pressure: $\pm 3$ mmHg<br>Pulse: $\pm 5\%$ of reading                                      |
| Measuring resolution:  | 1 mmHg                                                                                     |
| Inflation:             | Automatic inflation by internal pump                                                       |
| Memory function:       | 120 memory $\times$ 2 users (SYS, DIA, Pulse)                                              |
| Decompression:         | Constant exhaust valve system                                                              |

|                                  |                                                                                                      |
|----------------------------------|------------------------------------------------------------------------------------------------------|
| Power source:                    | 2- size »AAA« alkaline Batteries                                                                     |
| Rated voltage:                   | DC 3V 1.5W (direct current)                                                                          |
| Operation temperature:           | 5~40°C/41~104°F                                                                                      |
| Operation humidity:              | 15%~85%RH maximum                                                                                    |
| Storage temperature:             | -10~55°C/14~131°F                                                                                    |
| Cuff size:                       | 13.5-21cm                                                                                            |
| Storage humidity:                | 10%~95%RH maximum                                                                                    |
| Dimensions:                      | 70x65x30±1.0 mm                                                                                      |
| Device weight:                   | 129 g±5g (including batteries and cuff)                                                              |
| Batteries weight:                | 23 g±1.0g                                                                                            |
| Storage box weight:              | 60 g±1.0g                                                                                            |
| Size of the Sys /<br>Dia digits: | 12 mm                                                                                                |
| Time display:                    | 24 h display                                                                                         |
| Visible display size:            | 43.8x34 mm                                                                                           |
| Cuff pressure display<br>range:  | 0~290 mmHg                                                                                           |
| Electrical shock<br>protection:  | Internal power unit                                                                                  |
| Safety classifications:          |  Type BF equipment |

|                                      |                                                     |
|--------------------------------------|-----------------------------------------------------|
| Mode of operation:                   | Continuous operation                                |
| Protection against ingress of water: | IP22                                                |
| Accessories:                         | Storage case, 2 »AAA« batteries, instruction manual |

## 13. Manufacturer's Declaration

The Wrist Digital BP Monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the Wrist Digital BP Monitor should assure that it is used in such an environment. Electromagnetic Emissions: (IEC60601-1-2)

| Emission Test                                  | Compliance     | Electromagnetic Environment                                                                                                                                                                                                               |
|------------------------------------------------|----------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| RF emission<br>CISPR 11                        | Group 1        | The Wrist Digital BP Monitor uses RF energy only for internal functions. Therefore, this RF emission is extremely weak and there is little chance of it creating any kind of interference whatsoever with nearby electronic equipment.    |
| RF emissions CISPR 11                          | Class B        | The Wrist Digital BP Monitor is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes. |
| Harmonic emissions IEC 61000-3-2               | Not applicable |                                                                                                                                                                                                                                           |
| Voltage fluctuations/<br>flicker IEC 61000-3-3 | Not applicable |                                                                                                                                                                                                                                           |

### Electromagnetic Immunity: (IEC60601-1-2)

| Immunity test                                  | IEC60601-1-2 test level    | Compliance level           | Electromagnetic environment - guidance                                                                                                         |
|------------------------------------------------|----------------------------|----------------------------|------------------------------------------------------------------------------------------------------------------------------------------------|
| Electrostatic discharge (ESD)<br>IEC 61000-4-2 | ±6 kV contact<br>±8 kV air | ±6 kV contact<br>±8 kV air | Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %. |

|                                                                                                     |                                                                                                                                                                                            |                |                                                                                                                                                                                                                                                                                                         |
|-----------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Electric fast transient/ burst IEC 61000-4-4                                                        | $\pm 2$ kV for power supply lines<br>$\pm 1$ kV for input/output lines                                                                                                                     | Not applicable | Mains power quality should be that of a typical commercial or hospital environment.                                                                                                                                                                                                                     |
| Surge IEC 61000-4-5                                                                                 | $\pm 1$ kV differential mode<br>$\pm 2$ kV common mode                                                                                                                                     | Not applicable | Mains power quality should be that of a typical commercial or hospital environment.                                                                                                                                                                                                                     |
| Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11 | $< 5\% U_T$ (95% dip in $U_T$ ) for 0.5 cycle<br>$40\% U_T$ (60% dip in $U_T$ ) for 5 cycles<br>$70\% U_T$ (30% dip in $U_T$ ) for 25 cycles<br>$< 5\% U_T$ (95% dip in $U_T$ ) for 5 sec. | Not applicable | Mains power quality should be that of a typical commercial or hospital environment. If the user of the upper arm stlye requires continued operation during power mains interruptions, it is recommended that the Wrist Digital BP Monitor be powered from an uninterruptible power supply or a battery. |
| Power frequency (50/ 60 Hz) magnetic field IEC 61000-4-8                                            | 3 A/m                                                                                                                                                                                      | Not applicable | Not applicable                                                                                                                                                                                                                                                                                          |

Note:  $U_T$  is the a.c. mains voltage prior to application of the test level.

Note1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

- a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Wrist Digital BP Monitor is used exceeds the applicable RF compliance level above, the Wrist Digital BP Monitor should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the Wrist Digital BP Monitor.
- b) Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

| Immunity test                 | IEC60601-1-2 test level               | IEC60601-1-2 test level | Electromagnetic environment - guidance                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |
|-------------------------------|---------------------------------------|-------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Conducted RF<br>IEC 61000-4-6 | 3 Vrms 150 kHz to 80 MHz 80% AM (2Hz) | Not applicable          | Portable and mobile RF communications equipment should be used no closer to any part of the Wrist Digital BP Monitor, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.<br><br><b>Recommend separation distance</b><br>3V<br>$d = 1.2 \times P^{1/2}$ 80MHz to 800 MHz<br>$d = 2.3 \times P^{1/2}$ MHz to 2.5 GHz<br>where $P$ is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and $d$ is the recommended separation distance in meters (m).<br>Field strengths from fixed RF transmitters as determined by an electromagnetic site survey <sup>a</sup> , should be less than the compliance level in each frequency range <sup>b</sup> . Interference may occur in the vicinity of equipment marked with the following symbol:<br> |
| Radiated RF<br>IEC 61000-4-3  | 3 Vrms 80 MHz to 2.5 GHz 80% AM (2Hz) | 3 V/m                   |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |

## Recommended Separation Distances:

Recommended separation distance between portable and mobile RF communications equipment and the Wrist Digital BP Monitor.

The Wrist Digital BP Monitor is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Wrist Digital BP Monitor can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Wrist Digital BP Monitor as recommended below, according to the maximum output power of the communications equipment.

| Rated maximum output power of transmitter (W) | Separation distance according to frequency of transmitter m |                                               |                                                |
|-----------------------------------------------|-------------------------------------------------------------|-----------------------------------------------|------------------------------------------------|
|                                               | 150 kHz to 80 MHz<br>$d = 1.2 \times P^{1/2}$               | 80 MHz to 800 MHz<br>$d = 1.2 \times P^{1/2}$ | 800 MHz to 2.5 GHz<br>$d = 2.3 \times P^{1/2}$ |
| 0.01                                          | 0.12                                                        | 0.12                                          | 0.23                                           |
| 0.1                                           | 0.38                                                        | 0.38                                          | 0.73                                           |
| 1                                             | 1.2                                                         | 1.2                                           | 2.3                                            |

| Rated maximum output power of transmitter (W) | Separation distance according to frequency of transmitter m |                                               |                                                |
|-----------------------------------------------|-------------------------------------------------------------|-----------------------------------------------|------------------------------------------------|
|                                               | 150 kHz to 80 MHz<br>$d = 1.2 \times p^{1/2}$               | 80 MHz to 800 MHz<br>$d = 1.2 \times p^{1/2}$ | 800 MHz to 2.5 GHz<br>$d = 2.3 \times p^{1/2}$ |
| 10                                            | 3.8                                                         | 3.8                                           | 7.3                                            |
| 100                                           | 12                                                          | 12                                            | 23                                             |

For transmitters rated at a maximum output power not listed above, the recommended separation distance  $d$  in metres (m) can be determined using the equation applicable to the frequency of the transmitter, where  $P$  is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies

Note2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

### Remark:

|                                                                                   |                                                                                                          |                                                                                   |                                                                                                                                                           |
|-----------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------|
|   | Electrical devices are recyclable material and should not be disposed of with household waste after use! |  | Medical device                                                                                                                                            |
|   | Manufacturer                                                                                             |  | Unique Device Identifier                                                                                                                                  |
|   | Inapplicable for baby                                                                                    |  | Number of products in one packaging                                                                                                                       |
|   | Warning! Consult accompanying documents                                                                  |  | Before use, read the instructions. Electronic instructions for use: <a href="http://www.mediblink.com/f/m505.pdf">http://www.mediblink.com/f/m505.pdf</a> |
|   | Product reference number                                                                                 |  | Type BF equipment                                                                                                                                         |
|   | Model number                                                                                             | <b>IP22</b>                                                                       | Degrees of protection                                                                                                                                     |
|   | Batch number*                                                                                            |  | Authorised Representative in EU                                                                                                                           |
|  | Serial number                                                                                            |                                                                                   |                                                                                                                                                           |

\*The production date can be read from LOT number [YYYYMM]; first four digits represent the Year and last two digits Month of production. Example: LOT 202503 = March 2025

## 14. Warranty

**Product: Mediblink Wrist Blood Pressure Monitor M505**

**Manufactured for (Importer for EU & Distributor):** Mediblink d.o.o.,  
Gubčeva cesta 19, 8210 Trebnje, Slovenia; info@mediblink.com; www.mediblink.com

**Sellers name, address, signature and stamp\*:**

**Date of extradition/sales\*:**

\*If the invoice is accompanied by this warranty, and if all above information can be seen from the invoice, it is not necessary to fill in this field.

### WARRANTY TERMS

Dear customers!

The warranty period is **5 years** and starts from the day of product purchase. The warranty on cuff is **5 years**. In case of product claim, you have to show the invoice. We kindly ask you to save the invoice!

Unfortunately, wrong handling with the device is a reason for 95% of customer complains. You can easily avoid any problem, by getting useful information provided by our special service department. To reach our service department, you can call or send e-mail to Mediblink local distributor.

Before sending the product back to retailer, we kindly ask you to call our service department, to get help about how to use the device to save you with unneeded trips.

The manufacturer guarantees free elimination of all imperfections due to defects in material or manufacturing procedure by repairing or replacing the product. In case that the product can not be repaired or replaced, the customer will get the money refund. The guarantee is not valid in case of the force majeure, accidents or unexpected events (such as lightning, water, fire etc.), incorrect use or incorrect transport, non-compliance with safety and maintaining regulations or in case of unprofessional product intervention.

Traces of every day product usage (scratches, abrasions) and not subject to claim. The warranty does not eliminate the customer rights, which originate from seller responsibility for product flaws. By accepting the claimed product by the service department, the service department does not take responsibility for loss of saved data or settings on the product. All product repairs, which are performed out of product warranty period, have to be paid by customer by prior notice.

The manufacturer guarantees the product quality and flawless product operation in the warranty period, which starts with the day of product purchase. If the product can not be repaired in 45 days, the product will be replaced with a new one. In case that the product can not be replaced, the money will be refunded to the customer.

In case of product claim, call or send e-mail to Mediblink local distributor.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the Competent Authority of the Member State in which the user and/or patient is established.

**Replacement parts can be ordered at [www.mediblink.com](http://www.mediblink.com).**

**Manufacturer Information:**  
Shenzhen Combei Technology Co., Ltd.  
11-5B, No.105, Huanguan South Road,  
Dahe Community, Guanlan,  
Longhua District, ShenZhen,  
Guangdong, China

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